




## 100kN Wrap Grips

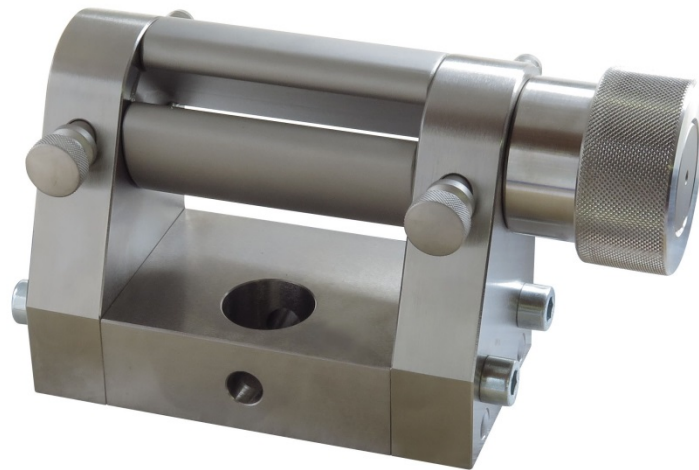

<u>Max Capacity:</u>	100kN/22K lbs.
<u>Roll Diameter:</u>	60mm
<u>Max Sample Width:</u>	80mm
<u>Max Sample Thickness:</u>	11mm
<u>Min Sample Length:</u>	910mm
<u>Weight:</u>	7.3kg each – steel.
<u>Coupling:</u>	Female 40mm w/ 16mm pin hole

Notes: The ribbon sample is put through the center of the drum and then wrapped by turning the knob. After a few turns of the knob two locking pins are used to secure the drum into place. Please take note to orient the grips so that the tensile force is translated through the center (See Below - Top Diagram).

## Schematic (mm)

